



- *Customized pressure ulcer prevention*
- *Patient safety and more dignity*
- *Improved working environment*

Careturner® - the unique solution for automatic pressure ulcer prevention and repositioning of patients



WHAT IS CARETURNER?

Careturner® is an innovative solution designed for full-automatic pressure relief and support of care staff with moving and handling of patients.

Careturner® can be operated by only one caregiver, when assisting with heavy lifts. It can be used for positioning or enveloping the patient to stimulate the senses. Patients can enjoy peaceful and undisturbed sleep during the customized pressure ulcer prevention program.

FEATURES:



Careturner works in conjunction with the nursing bed to optimize comfort.



Careturner can cradle and "envelop" the patient giving comfort and stimulating senses, when both wings are elevated.



The wings can be elevated up to 80° for moving and handling purposes, such as personal care and sling application




Programmable control box

Automatic use of Careturner®

The programmable control box has an option for two customized programs chosen by the caregiver, ensuring optimal pressure ulcer prevention for the patient. Careturner is turning the end-user gently from side to side during the night giving pressure relief and sense stimulation. Careturner operates very silent and without any vibration allowing the patient to sleep undisturbed. In automatic mode the wings can elevate from 0-40°.

Manual use of Careturner®

Operated by the care staff. With a push on a button it turns the end-user from one side to the other to support the care staff in daily moving and handling procedures e.g. sling application and personal care. In manual mode the wings can elevate from 0-80°.

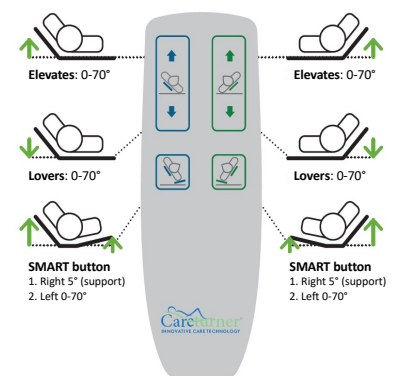

Hand control

CARETURNER RELEASES RESOURCES IN THE CARE SECTOR



3-IN-1 SOLUTION BENEFITS WITH CARETURNER®



BETTER COMFORT AND UNDISTURBED SLEEP

- Prevents pressure ulcers
- More peace of mind and less pain
- Better sleep
- Fewer conflicts between patient and caregiver



BETTER WORKING ENVIRONMENT

- Safe patient handling
- Fewer heavy lifts
- Reduced attrition of care staff
- Fewer work related sick days



BETTER RESSOURCE UTILIZATION


- Automated turning saves resources
- Efficient mobilization (from 2 to 1 caregiver)
- Time savings for the organization (ROI: 3-4 months)
- Smarter processes free up resources for better quality in care

13.6 million
shortfall of formally
employed
caregivers world wide
by 2030

** World Health Organization*

Forecast for the development of care
needs worldwide*

Number of people in millions

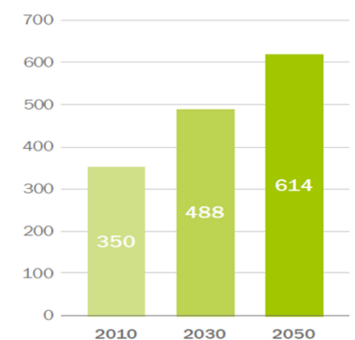

THE COST OF PREVENTION OF PRESSURE ULCERS IS FAR LOWER THAN IN TREATING PRESSURE ULCERS

4 million cases of people suffering from pressure ulcers in Europe*

It is estimated that 95% of all pressure ulcers are preventable

Prevention of pressure ulcers costs per patient pr. Day up to 87€

Treatment of pressure ulcers costs per patient pr. day up to 470€

**Source: EPUAP & National Clinical Guideline for the Prevention of Pressure Ulcers*



Prevention of pressure ulcers with Careturner® for use in the fields of

- Intensive care
- Care for residents affected by dementia
- Palliative care
- Post-acute care and rehabilitation
- Geriatric psychiatry
- Preventive care
- Bariatric patients

Careturner Case Studies



Careturner®
INNOVATIVE CARE TECHNOLOGY

Click

Case Studies

Careturner A/S - Lyskær 8B
2730 Herlev - Denmark

+45 20978115
info@careturner.com
www.careturner.com



CARETURNER® CAN BE USED IN COMBINATION WITH

 **wissner-**
 **bosserhoff**

SENTIDA 5



SENTIDA 6



SENTIDA 6 HOSPITAL



Angles in manual mode:

0-80°



Angles in automatic mode:

0-40°



Time interval in automatic mode:

0-180 min



Max user weight:

**PW 165 kg
SWL 185 kg**



Heaviest part of the product:

18 kg



Degrees the blue cover can be washed:

80°