

The unique solution for automatic and personalized pressure ulcer prevention, safe patient handling and sensory stimulation



- Personalized pressure ulcer prevention
- Sensory stimulation enhancing sleep and well-being
- Improved working environment and safe patient handling
- Frees up time for better care

What is Careturner?

Careturner® is an innovative solution designed for full-automatic personalized pressure ulcer prevention, sensory stimulation and safe patient handling.

FEATURES:



Personalized Pressure Ulcer Prevention

Careturner is controlled by a computer and does not make any noise or vibration, ensuring comfortable and gently relief of pressure points when repositioning the patient from side to side.



Sensory Stimulation

Careturner can also elevate both wings, enveloping the patient which has a calming effect and optimizing comfort, well-being and improving the sleep.



Safe Patient Handling

The Careturner can be operated manually, making it ideal for moving and handling within single carer working environments.



Follows The Beds Movements

Careturner is easy to install and ensures that all the functions of the bed can still be used for optimal comfort.



Automatic use of Careturner®


Personalize two pressure ulcer prevention programs. Time and angles for each side can be defined. Careturner operates very silent and without any vibration, allowing the patient to sleep undisturbed.

- Time from 0-180 min. between each change of position
- Angles from 0-40°
- Possible to relief of one side if necessary (left, right or back)
- Counterhold from 0-12° the opposite side reduces shear and improves comfort
- Timer function allows Careturner to start at desired time
- Switch from automatic to care function without waiting time
- 30 min. pause function allowing manual usage (switches back to a automatic use)

Manual use of Careturner®

Ideal use for safe patient handling, when turning a patient manually, applying a sling or by personal care.

- Operate the wings simultaneously or alone
- Wings can be elevated from 0-80°
- Smart buttons ensure automatic counterhold from the opposite side to reduce shear
- Reduces unwanted physical contact with the patient, ideal for extraverted patients
- Smart buttons move both wings simultaneously to ensure maximum comfort and safety

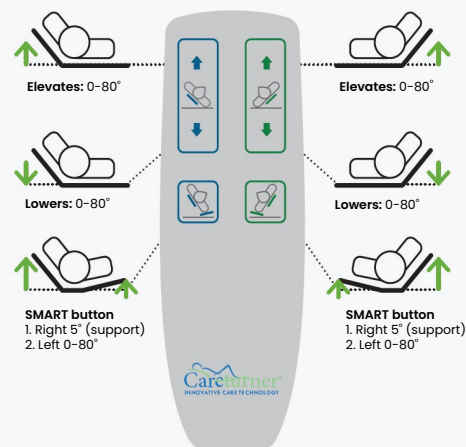

Treatment of pressure ulcers is 5 times more expensive than prevention

Pressure mapping 30° with 9° counterhold



Relief of right side



Careturner not activated



Relief of left side

Pressure mapping "the hug"




Careturner not activated



15-15°



20-20°

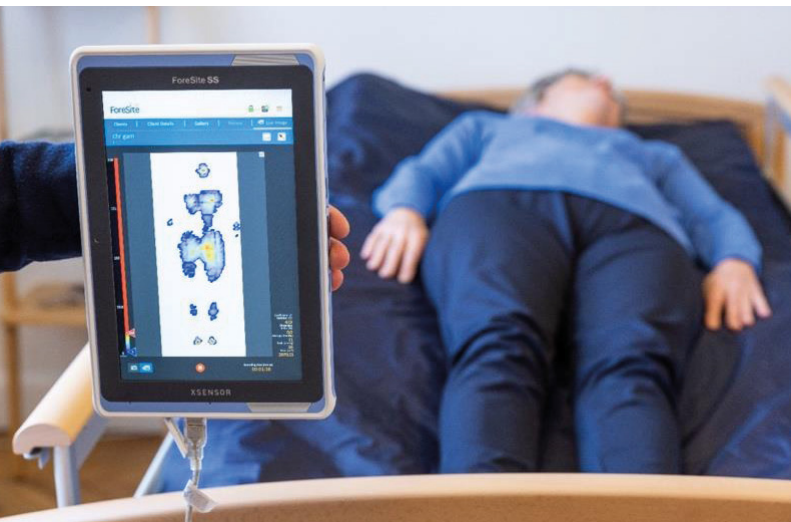

Latest Health Technology Assessments documents:

95% of all patients should not be repositioned anymore after using The Careturner.

50% of all patients feel more calm in general and are easier to work with when moving and handling in bed.

The Careturner is considered a technology which can attract new and keep current healthcare workers.

4 million cases of people suffering from pressure ulcers in Europe*

It is estimated that 95% of all pressure ulcers are preventable
Prevention of pressure ulcers costs per patient pr. day up to 87€
Treatment of pressure ulcers costs per patient pr. day up to 470€

**Source: EPUAP & National Clinical Guideline for the Prevention of Pressure Ulcers*

Prevention of pressure ulcers with Careturner® for use in the fields of



- Dementia
- Home care
- Palliative care
- Acute care
- Rehabilitation
- Geriatric psychiatry
- Preventive care
- Intensive care

Careturner releases resources in the healthcare sector

3-IN-1 SOLUTION BENEFITS WITH CARETURNER®



BETTER COMFORT AND UNDISTURBED SLEEP

- Prevents pressure ulcers
- More peace of mind and less pain
- Better sleep and improved comfort
- Fewer conflicts between patient and caregiver



IMPROVED WORKING ENVIRONMENT

- Safe patient handling
- Fewer heavy lifts
- Reduced attrition of care staff
- Fewer work related sick days



BETTER RESOURCE UTILIZATION


- Automated turning saves resources
- Efficient mobilization (from 2 to 1 caregiver)
- Time savings for the organization (ROI: 3-4 months)
- Smarter processes free up resources for better quality in care
- Reduces the number of outings for home care

13.6 million
shortfall of
caregivers world
wide by 2030

* World Health Organization

Forecast for the development of care needs worldwide*

Number of people in millions

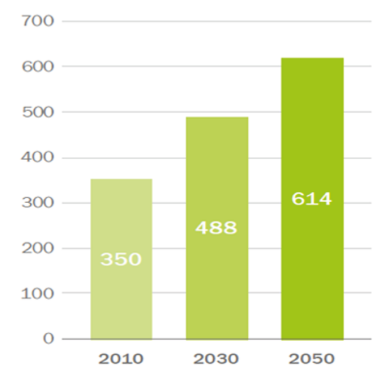

CARETURNER® CAN BE USED IN COMBINATION WITH:



Angles in manual mode:

0-80



Angles in automatic mode:

0-40



Time interval in automatic mode:

0-180 min



Max user weight:

**PW 165 Kg.*
SWL 185 Kg.***



Heaviest part of the product:

18 kg



Degrees the blue cover can be washed:

80

*Specifications depend on the individual models



Careturner Academy

Scan for tutorials, guides and documentation



Careturner A/S - Lyskær 8B, 2730 Herlev - Denmark +45 20978115 info@careturner.com www.careturner.com